beurer

BC

SL 70

List of supported
devices

READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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CONTENT

1. IMPORTANT SAFETY INSTRUCTIONS.....	3
2. Introduction	4
3. Parts and Controls.....	6
4. To Set Up	6
5. To Use	7
6. Care and Maintenance	8
7. Specifications	8
8. FCC Compliance	9
9. Warranty	9

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

The following signs appear in the Safety Section (page 3).

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.

CAUTION

CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

This product satisfies the requirements of the applicable European and national directives.

Disposal in accordance with EC Directive WEEE (Waste Electrical and Electronic Equipment).

CAUTION – To reduce the risk of personal injury or product/property damage:

- The device is intended solely for the purpose described in these operating instructions. The manufacturer is not liable for damage occurring as a result of inappropriate or careless use.

- Protect the device from impacts, humidity, dirt, marked temperature fluctuations and direct sunlight.
- Keep packaging material and the device away from children.
- Keep the device away from heat sources and do not use near shortwave or microwave devices.
- Use and store the snore stopper at an ambient temperature between 32° F and 140° F (0° C - 60° C).
- The snore stopper is splash-resistant only. Keep the device away from water.
- Do not use any chemical or abrasive cleaning agents. Clean the device occasionally with a 70% isopropyl alcohol pad.
- Fully charge the battery at least every three months.

Rechargeable Battery Handling Safety Precautions

- If your skin or eyes come into contact with fluid from the battery cell, flush out the affected areas with water and seek medical assistance.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect rechargeable batteries from excessive heat.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush rechargeable batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter “4. Charging the Battery”).
- In order to achieve the longest battery service life possible, fully charge the battery at least every three months.


2. Introduction

A restorative night's sleep is very important for our physical and psychological well-being; during sleep we process things such as new experiences and improve our mental fitness. If you snore, however, your

sleep quality is affected. This device is designed to reduce snoring and increase sleep quality as well as your well-being, so you can get a better night's sleep.

How does the snore stopper work?

When attached to the left or right ear, the snore stopper detects snoring noises, which produce vibrations transmitted via the cranial bones. When the snore stopper detects noises and vibrations, it emits impulses to the ear that consist of tones and/or vibrations. The intensity

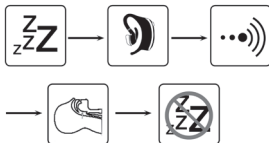of these impulses starts at a low level and increases until snoring stops. The throat muscles become firm and the airways open as a reaction to these impulses, therefore calming breathing and decreasing snoring. In the long term, this training of the muscles in the throat produces a lasting reduction in snoring. Note that it may take up to seven days to become accustomed to the device.

Download and install the free 'beurer SleepQuiet' app to access the following functions:

- Detailed snoring analysis including snores per hour, recording duration, and maximum snoring volume
- Calendar with snoring diary
- Evaluation view (in days/weeks/year)
- Adjustable intensity and sensitivity

App Requirements:

- iOS version 8.0 or later, or
- Android™ version 4.4 or later
- Bluetooth® 4.0

For a list of supported devices, scan this QR code:

3. Parts and Controls

- 1 ON/OFF button
- 2 Operation light
Green = battery at 41–100%
Red = battery at 0–40%
- 3 *Bluetooth*[®] light
- 4 Micro USB connection
- 5 Earplugs

4. To Set Up

Using with the ‘beurer SleepQuiet’ app

Using this device with the ‘beurer SleepQuiet’ app is not required, but the app does provide additional features such as data recording, statistics, and a calendar function, among others. To pair your mobile device with the Snore Stopper:

1. Ensure that your smartphone is connected to the internet.
2. Activate *Bluetooth*[®] on your mobile device.
3. Download and install the free ‘beurer SleepQuiet’ app from either the Apple App Store (for iOS) or Google Play (for Android devices).

-  Depending on the settings for your smartphone, additional charges may occur for connecting to the Internet or transferring data.
4. Start the ‘beurer SleepQuiet’ app and follow the instructions.

Charging the Battery


1. Lift the cover from the micro USB port.
2. Connect the included charging cable to the micro USB port and to either an available USB port on your PC (ensure PC is powered on) or

to the included AC Adapter. If using the AC Adapter, plug the adapter into a 110-VAC electrical outlet.

3. While charging, the operation light flashes green and will illuminate steadily in green when the battery is fully charged.
4. The battery should be charged in two hours. At that time, disconnect the charging cable from the computer or AC Adapter and the micro USB port.
5. The device will operate 4-5 days on a full charge. When the operation light turns red, the battery level is at 41% or less.

5. To Use

1. Press and hold the ON/OFF button for three seconds; the operation light will illuminate. The LED will shut off after a short time so it does not disturb you while sleeping.
2. Place the device onto either ear; the right ear is shown here. Do not bend the upper clip too much to avoid damaging the device. Turn the earplug so that it faces into the ear canal as shown. If the earplug moves around easily, try a larger one.

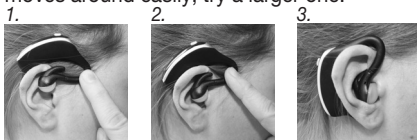3. Go to sleep. When you wake up, shut the device off by pressing and holding the ON/OFF button for about three seconds; the device will vibrate briefly to indicate it is switched off.

Transferring data to your mobile device

To transfer the recorded data from your snore stopper to the 'beurer SleepQuiet' app, proceed as follows:

1. Activate *Bluetooth®* on your mobile device.
2. Open the 'beurer SleepQuiet' app. Press the ON/OFF button on the snore stopper; recorded data will transfer to the smartphone. A large volume of data may take some time.

6. Care and Maintenance

- Shut the device off before cleaning.
- Clean the device with a soft damp cloth then wipe dry with a soft, lint-free cloth after each use.
- Only clean the device as instructed. Always keep liquids from entering the device or accessories.
- Do not use the device until it is completely dry.
- An elongated object or brush works well to remove any ear wax stuck in the earplug cover.
- Allow the device to dry fully after cleaning. Store the device in a safe, dry and clean place.

7. Specifications

Dimensions:	1.9 in x 0.6 in x 2.4 in (49 x 15 x 62 mm)
Weight:	15 g (without earplugs)
Connection:	Micro USB
Power supply:	Input: 100–240 VAC; 50/60 Hz; 0.3 A Output: 5 V $\overline{\text{---}}$ 1000 mA
Battery:	80 mAh lithium polymer
Battery Usage:	About four days on a full charge
Memory:	One night of heavy snoring
Earplug material:	Silicone
Data transfer:	<i>Bluetooth</i> ® low energy technology
Operating conditions:	Temperature: 32°–140° F (0°–60° C) Relative humidity: 30%–75% Ambient pressure: 860–1060 hPA
Storage and transportation conditions:	Temperature: -4°–140° F (-20°–60° C) Relative humidity: < 75% Ambient pressure: 860–1060 hPA

The *Bluetooth*® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer GmbH is under license. Other trademarks and trade names are those of their respective owners.

8. FCC Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

9. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Snore Stopper, Model SL70, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Snore Stopper, Model SL70, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement. For warranty service contact our customer service department at **1-800-536-0366 or at info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/

or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

Distributed by:
Beurer North America LP,
900 N Federal Hwy Ste 300
Hallandale Beach, FL 33009 USA
Made in China

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The *Bluetooth®* word mark and logos are registered trademarks owned by Bluetooth SIG. Inc. and any use of such marks by Beurer GmbH is under license. Other trademarks and trade names are those of their respective owners.